











Simple fitting supportive sitting with built in pressure care*



Hours sitting in a poorly fitted chair is hard work! The Avanti can easily be adjusted to provide the necessary postural support to maximise comfort and relaxation.

Correct leg length, arm height, seat dimensions, back and seat angles provides stability, reduces point pressures and counters slumping over

*LowzoneTM Contour cushion as standard

Our contoured LowzoneTM viscoelastic cushions are regarded as providing "Medium High" pressure relief.

In our own trials using our FSA pressure monitors, these regularly exceed non contoured gel, water and standard viscoelastic cushions in pressure reduction. They are inherently more stable than fluid based systems.

Tilt in space

Forward for function or access Back to offload the bony prominences and to relax

Independently adjustable back

Allows you to get the final fit right.

Independently adjustable leg rest

For a more relaxed position



Optional Head cushion



Quick and Easy Footplate adjustment

Options

- Extra lateral support
- Head cushion
- Footplate for non ambulant users
- Tray









<u>Adjustments</u>

| Seat Depth | | 380-550mm |
|--------------------------------|------------------|----------------|
| Seat Width - Snug arm | | 380-450mm |
| Seat Width - Std. arm | | 450-550mm |
| Back height Std. | | 830mm |
| Back Height Tall | | 900mm |
| Seat to Floor (Std. no footpl | ate) | 337-437mm |
| Seat to Floor (Tall no footpla | ate) | 337-437mm |
| Seat to Footplate - Std. | | 257-357mm |
| Seat to Floor - Tall | | 382-482mm |
| Tilt in Space - | Manual Operation | |
| Back recline - | Manual Operation | |
| Legrest - | Manual Operation | |
| Footplate - | Mar | nual Operation |
| | | |

Fabric: Cuero

Soft, vapour permeable waterproof PU


Black

Blue

Green

Red

Brown

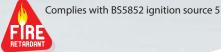







The Helping Hand Company Bromyard Road, Ledbury Herefordshire HR8 1NS Tel: **01531 635388**

www.**Symmetrikit.com** sales@Symmetrikit.com

